MAGUS CHD30 DIGITAL CAMERA

MAGUS

www.magusmicro.com

Levenhuk Inc. (USA):
928 E 124th Ave. Ste D, Tampa, FL 33612, USA,
+1-813-468-3001,
contact_us@levenhuk.com


Levenhuk Optics s.r.o. (Europe):
V Chotejně 700/7, 102 00 Prague 102, Czech Republic,
+420 737-004-919,
sales-info@levenhuk.cz

MAGUS® is a registered trademark of Levenhuk, Inc.

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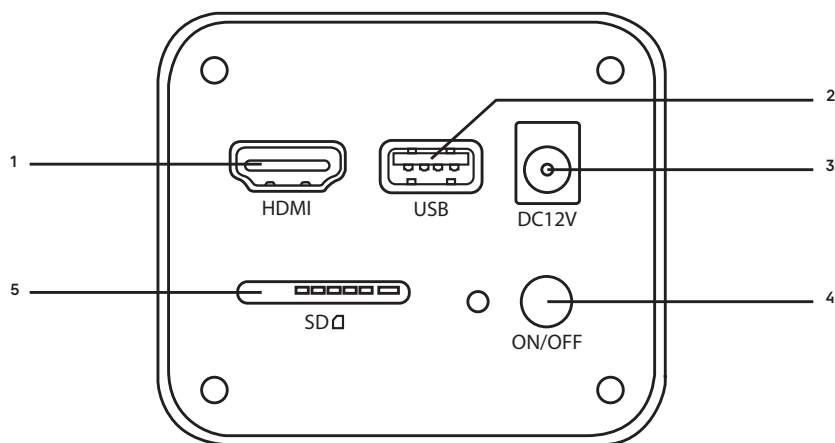EN MAGUS CHD30 DIGITAL CAMERA

Never look at the sun, a bright light source, or laser beam through the camera. It can cause permanent damage to your eyes.

MAGUS high-speed HDMI cameras are designed to be paired with a digital microscope. Sony Starvis high-sensitivity sensor ensures perfect color reproduction and image quality even in low light conditions. You can install the camera software to view, save, and edit captured images and videos. The rugged aluminum housing is designed for operation at low temperatures (-10°C).

Before starting the camera please connect the standard C-mount camera to the camera adapter and connect it to the third tube of the microscope which will convey the middle image of the object to the camera sensor.

CONNECT THE CAMERA TO A PC/LAPTOP VIA WI-FI

Make sure your PC supports Wi-Fi.

1. Install the MAGUSView software from the USB flash drive (included in the kit) on the PC.
2. Connect the camera to the power supply (3). Turn the camera on (4).
3. Connect the Wi-Fi antenna to the connector (2).
4. After the antenna starts blinking, select the Wi-Fi network that you want to connect to on your PC.
The password is 12345678.
5. Open the MAGUS View program. Choose a camera from the list.

HDMI


1. Use an HDMI cable to connect the camera to an HDMI monitor (1).
2. Connect the USB mouse to the connector (2).
3. Insert the memory card into the memory card slot (5).
4. Connect the camera to the power supply (3). Turn the camera on (4).
5. Turn the monitor on and view the video in the built-in program.

CALIBRATION USING A CALIBRATION SLIDE

You should use a special stage micrometer (calibration slide) to determine the linear dimensions (in millimeters or microns).

The calibration slide is a transparent glass (of the same size as the specimen slide) that has a micrometer scale with a scale division of 0.01mm etched on the surface.

The calibration slide (stage micrometer) is used to calibrate the image analysis software for measurements in actual units. In the calibration mode, you should capture an image of the micrometer scale with every objective magnification and indicate the known distance. That lets you establish a scale of the image in actual units (micrometer, millimeter, etc.).

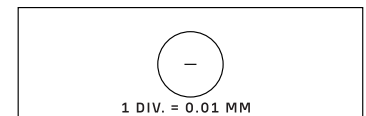Calibration:

1. Place the calibration slide on the microscope stage.
2. Select the required objective and set the maximum camera resolution.
3. Get a contrast image of the scale on the monitor screen and capture the image.
4. Select the "Calibrate" function in the software you are using.
5. Double-click on the maximum visible distance and enter the value in actual units.
6. Enter the calibration setting and check the result. The program will save the calibration factor.
7. You can select any measurement unit later, and all the results will be re-calculated in accordance with this selection.

CAMERA CONTROL PANEL

To open the Camera control panel, move the mouse cursor to the left of the window. The control panel will automatically pop up.

CAMERA CONTROL PANEL	FUNCTION	FUNCTION DESCRIPTION
	Snap	Captures the image and saves it to the SD card.
	Record	Records a video and saves it to the SD card.
	Auto Exposure	Automatically adjusts the exposure time and gain according to the exposure compensation value.
	Exposure Compensation	Available when Auto Exposure is selected. Adjusts the exposure compensation according to the current video brightness to achieve the proper brightness value.
	Exposure Time	Available when Auto Exposure is not selected. Adjusts the time of exposure and thus adjusts the brightness of the video.
	Gain	Adjusts the brightness of the video. The Noise will be adjusted accordingly.
	Red	Adjusts the proportion of the red color in RGB in the video.
	Green	Adjusts the proportion of the green color in RGB in the video.
	Blue	Adjusts the proportion of the blue color in RGB in the video.
	Auto White Balance	Adjusts the white balance of the video continuously. White balance is adjusted to obtain a natural ratio of cold and warm colors on the image with regard to the light source type.
	Sharpness	Adjusts the Sharpness level of the video.
	Denoise	Removes the noise from the video.
	Saturation	Adjusts the Saturation level of the video.
	Gamma	Allows you to adjust the midtone values of the image. The higher the Gamma value is, the darker the image will be.
	Contrast	Adjusts the Contrast level of the video.
	DC (Direct current)	Allows you to use a light source without light flickering.
	AC(50HZ) (Alternative current)	Eliminates light flickering when a 50Hz light source is used.
	AC(60HZ) (Alternative current)	Eliminates light flickering when a 60Hz light source is used.
	Default	Restores the default settings.

MEASUREMENT TOOLBAR

To open the Measurement Toolbar, move the mouse cursor to the top of the window. The toolbar will automatically pop up.

ICON	FUNCTION	ICON	FUNCTION
	Float/Fix switch of the Measurement Toolbar		Ellipse
<input checked="" type="checkbox"/> Visible	Show/Hide Measurement Objects		Annulus
Pixel	Select the desired Measurement Unit		Two Circles and its Center Distance: drawing up of two circles and the measurement of the distance between the centers of these two circles.
NA	Choose the same magnification value as the magnification value of the microscope to ensure accuracy of measurement results when the measurement is not in pixels.		Arc
	Object Select		Polygon
	Point		Curve
	Angle		Executes Calibration with the help of a micrometer to determine the corresponding relation between magnification and resolution as well as between the measurement unit and the sensor pixel size.
	Arbitrary line		Conjugated Calibration: press the CC button to perform the conjugated calibration. It will return the position of the camera sensor to the initial C-mount position. Then manually adjust the coarse and fine focus of the microscope. Make sure that the magnification value of the camera software corresponds to the magnification value of the microscope. Then select the measurement unit to perform the measurement.
	Parallel lines		Exports the Measurement information to CSV file (*.csv)
	Horizontal Line		Deletes all of the measurement objects
	Vertical Line		Settings
	Rectangle		Exits from Measurement mode
	Circle		Use these icons to select Move Left, Move Right, Move Up, Move Down, Color Adjustment, or Delete.

Select a figure to draw up; click at the desired initial point on the image, click at the desired end point. The selected figure will be drawn up on the image and the figure parameters will be displayed.

CAMERA CONTROL TOOLBAR


To open the Camera Control Toolbar, move the mouse cursor to the bottom of the window. The control toolbar will automatically pop up.

ICON	FUNCTION	ICON	FUNCTION
	Zoom In the Video Window		Zoom Out the Video Window
	Horizontal Flip		Vertical Flip
	Video Freeze		Displays the cross line
	WDR		Open the Auto Focus Control Panel
	Browses the images and videos on the SD Card		Settings
	Check the Version of MAGUSView		


AUTOFOCUS CONTROL PANEL

Press the **AF** button to open the AutoFocus Control panel in the right part of the window and the rectangle Focus Region in the center of the image.

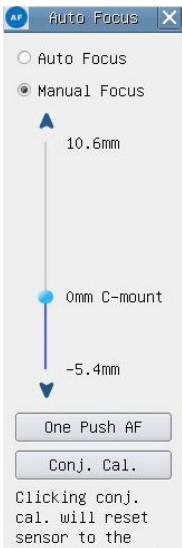Auto Focus	When the Auto Focus function is selected, the system will start autofocusing on the object.
Manual Focus	When the Manual Focus function is selected, you need to adjust the camera sensor position manually by using the mouse to scroll the mouse wheel up and down until the object is in focus.
One Push AF	Click the One Push AF button to perform the autofocusing operation just once.
Conjugated Calibration	Click the Conjugated Calibration button to reset the camera sensor to the standard C-mount position. The Conjugated Calibration function helps to calibrate the camera sensor position when the camera video display window is clear and the view from microscope eyepiece is clear. Apply the Conjugated Calibration function when starting the camera for the first time to set the camera sensor at the standard C-mount position. It is required to set the object plane, the eyepiece image plane and the camera adapter image plane at their standard positions. If the height of the object changes, you must ensure that the camera sensor is set at the standard C-mount position. Adjust the coarse and fine focus of the microscope to focus on the object. Prior to perform measurements, please start the Conjugated Calibration option to ensure the accuracy of the measurement results (see the Conjugated Calibration option in the Measurement Toolbar Section).

FOCUS REGION IN THE VIDEO WINDOW

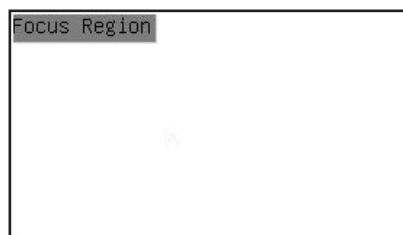
The Focus Region is used to select the region of interest (ROI) for the autofocus operation. You can click any part of the video display window to change the focus region position for the Auto Focus operation.

When you close the Auto Focus Control Panel, the Focus Region will also be closed automatically.

If the Focus Region function is enabled, the Measurement Toolbar will not appear when the mouse cursor is moved to top of the video display window.

SETTINGS

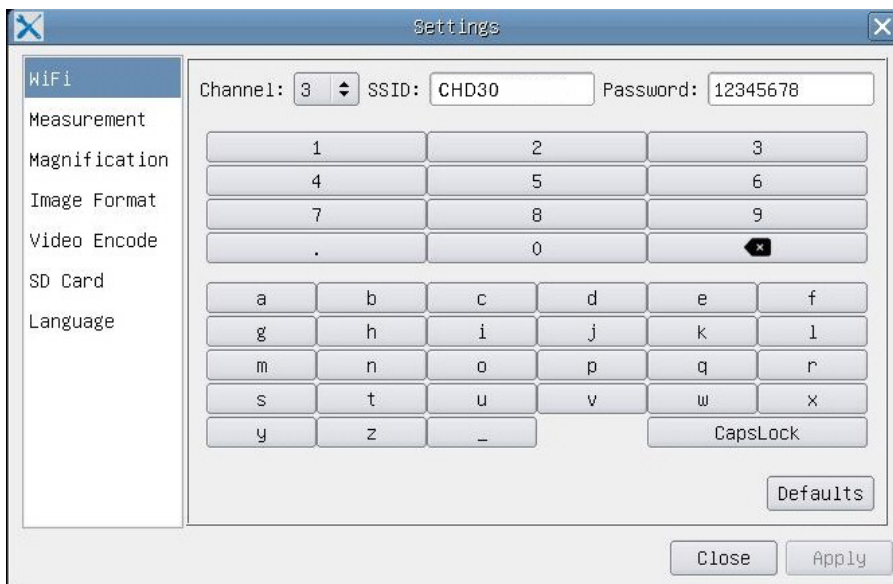
WI-FI

Channel: Avoids interference caused by using the same WiFi signal Channel. Suggests different channels for different cameras when several WiFi cameras are operated at the same time.

SSID: You can use the soft keyboard below to set the name of the WiFi signal.

Password: You can use the soft keyboard below to set the Password for the WiFi signal.

Default: Resets the Channel, SSID, Password options to the default values.

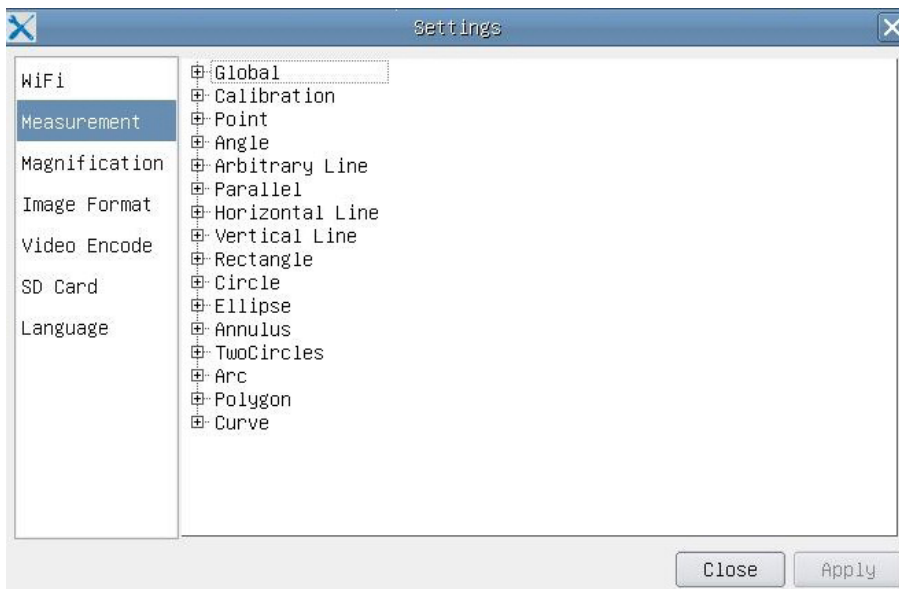MEASUREMENT

Global	Precision	Sets the required number of digits after the decimal point. The measurement result value will be displayed with the required precision.
Calibration	Line Width	Defines width of the lines for calibration.
	Color	Defines color of the lines for calibration.
	EndPoint	Type: defines shape of the endpoints of lines for calibration: Null means no endpoints, Rectangle means rectangle shape of endpoints. It facilitates the object alignment.
Point, Angle, Arbitrary Line, Parallel, Horizontal Line, Vertical Line, Rectangle, Circle, Ellipse, Annulus, Two Circles, Arc, Polygon, Curve		

You can Left-click the  and select Measurement to open the corresponding window of settings. You can specify the parameters of the required measurement objects.

MAGNIFICATION

Name	Sets the magnification item name when calibration, such as 4x, 10x, 100x, etc. Other information could be added too, for example, microscope model, operator name, etc.
Resolution	Sets the pixels per meter.
Clear All	Clears all of the calibrated magnification items.
Delete	Deletes the selected magnification item.

IMAGE FORMAT

Image Format	JPEG: you can save files in the JPEG format to get a high compression and high quality images and to save the storage space. TFT: you can save not only the images but also you can save the measurement data.
Measurement Object Saving Method	Burn-in Mode: merges the measurement objects into one current image. Further edition of the measurement objects is impossible. The changes are irreversible. Layered Mode: saves the measurement objects at different layers of the current image into the target file. The measurement objects in the target file can be edited with the use of corresponding software on the PC. The changes are reversible.

VIDEO

Fast Forward/Reverse Interval	Sets the time interval of the playback of video files.
Video Encode	You can select the video encoding format MJPEG or H264. MJPEG: Ensures high quality of image, but requires more storage space. H264: Ensures good quality of image and saves the storage space.

STORAGE

Storage Device	SD Card is the only supported storage device.
File System Format of the Storage Device	You can select the file system format: FAT32 or NTFS. Use PC to format the SD cards and switch between the file system formats. FAT32: The maximum video file size of a single file is 4GB. NTFS: The maximum video file size of a single file is 2048GB. Unknown Status: SD card not detected or the file system is not identified.

You can also select a language from the list.

THE KIT INCLUDES

MAGUS microscope digital camera, 12V, 1A power adapter, HDMI cable, USB mouse, 32GB SD card, LM816 USB Wi-Fi adapter, USB flash drive with drivers and software, mounting plate and screws, user manual, and a warranty card.

SYSTEM REQUIREMENTS

Windows 8/10/11 (32bit and 64bit), Mac OS X, Linux, up to 2.8GHz Intel Core 2 or higher, minimum 4GB RAM, USB2.0 port, 19" or larger display.

SPECIFICATIONS

	MAGUS CHD30
Maximum resolution	1920x1080px
Megapixels	2
Sensor	1/1.9" (7.20x4.05mm) Sony CMOS sensor
Pixel size	3.75x3.75µm
Light sensitivity	1120mv with 1/30s
Exposure	0.06–918ms
Video recording	+
Frame rate, fps at resolution	60@1920x1080px (HDMI) 25@1920x1080px (Wi-Fi)
Image format	*.jpeg
Video format	*.asf format with 1080P resolution

Shutter type	ERS
Auto focus	+
Interface	HDMI, Wi-Fi
USB Wi-Fi adapter	802.11n 150Mbps
Body	metal
Operating temperature range	-10...50°C

The manufacturer reserves the right to make any changes to the product lineup and specifications or discontinue the product without prior notice.

CARE AND MAINTENANCE

- Do not disassemble the camera on your own.
- Keep the camera away from moisture; do not use it in rainy weather.
- Keep the camera away from shock or excessive pressure.
- Do not overtighten the locking screws.
- Keep the camera away from hazardous environment, home and car heaters, incandescent lamps or open fire.
- When cleaning the lenses, first blow any dust or debris off the surface or wipe them off with a soft brush. Then wipe the lens with a soft cleaning cloth slightly damp with spirits or ether.
- Seek medical advice immediately if a small part or a battery is swallowed.

MAGUS WARRANTY

MAGUS optic products carry a **5-year warranty** against defects in materials and workmanship. All MAGUS accessories are warranted to be free of defects in materials and workmanship for **six months** from the purchase date. The warranty entitles you to the free repair or replacement of the MAGUS product in any country where a Levenhuk office is located if all the warranty conditions are met.

For further details, please visit: www.magusmicro.com

If warranty problems arise, or if you need assistance in using your product, contact the local Levenhuk branch.

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© 2006–2024 Levenhuk, Inc. All rights reserved.